



MICRO-Pedi™

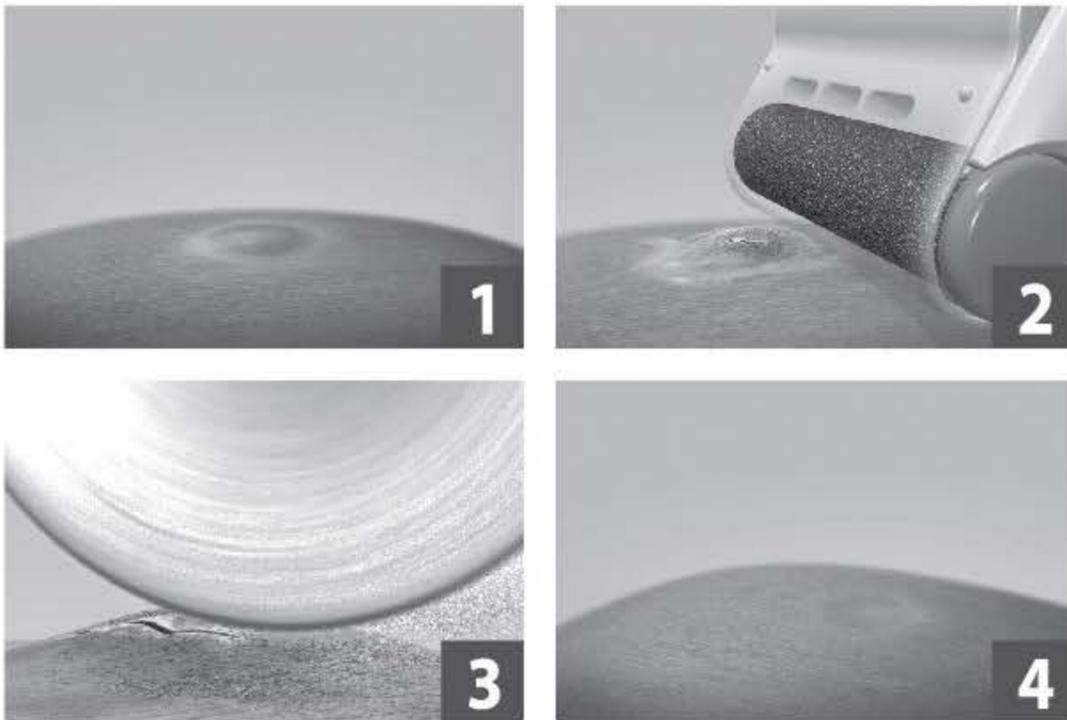


AP-3RP
Battery Operated

GENERAL DESCRIPTION OF PARTS

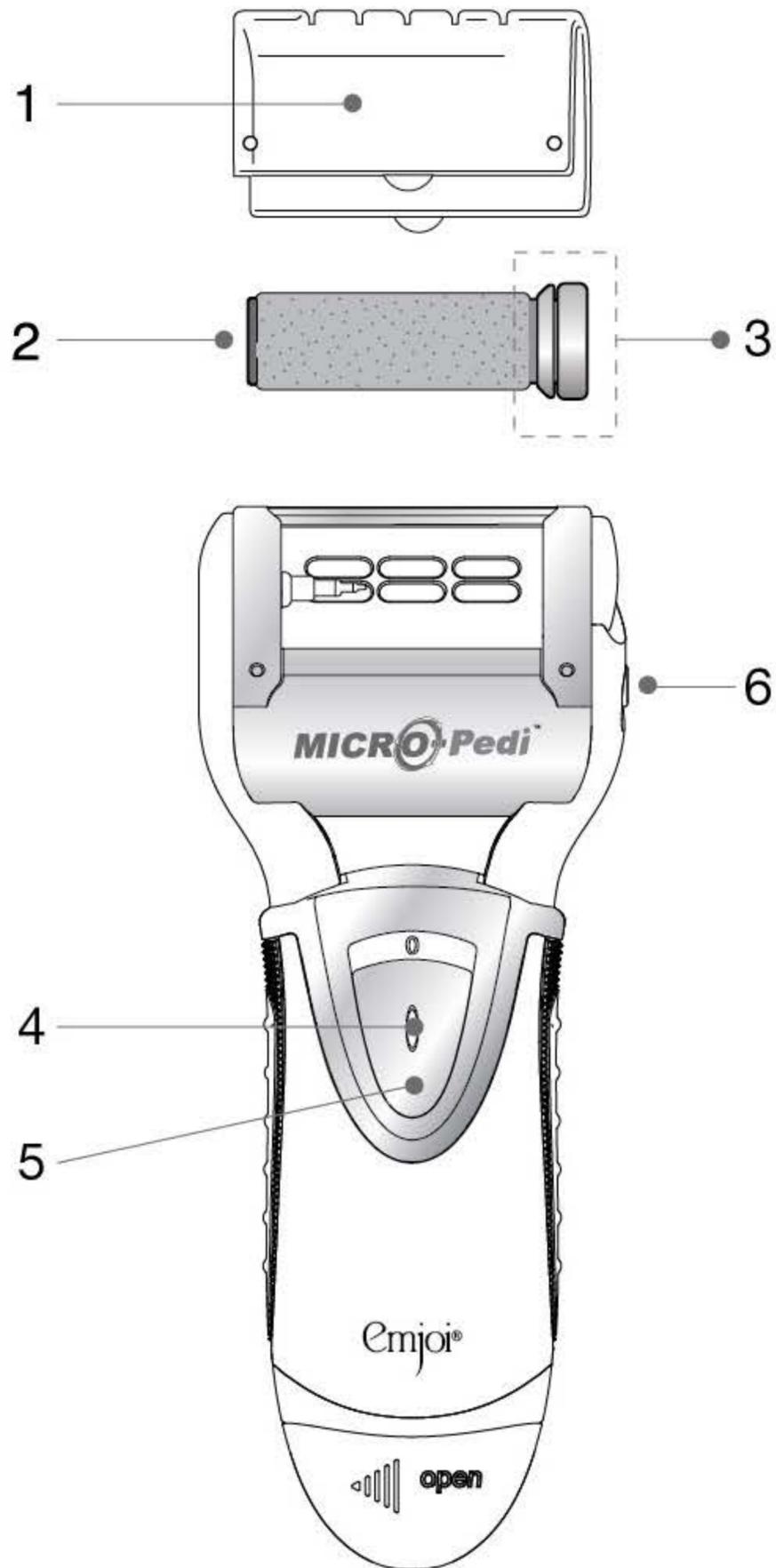

1. Protective Cover
2. Micro-Pedi Roller
 - blue for regular skin
3. End Cap of Micro-Pedi Roller
4. Safety Switch Lock
5. On/Off switch
6. Release Button for Micro-Pedi Roller
7. Battery cover (bottom cap of handle)

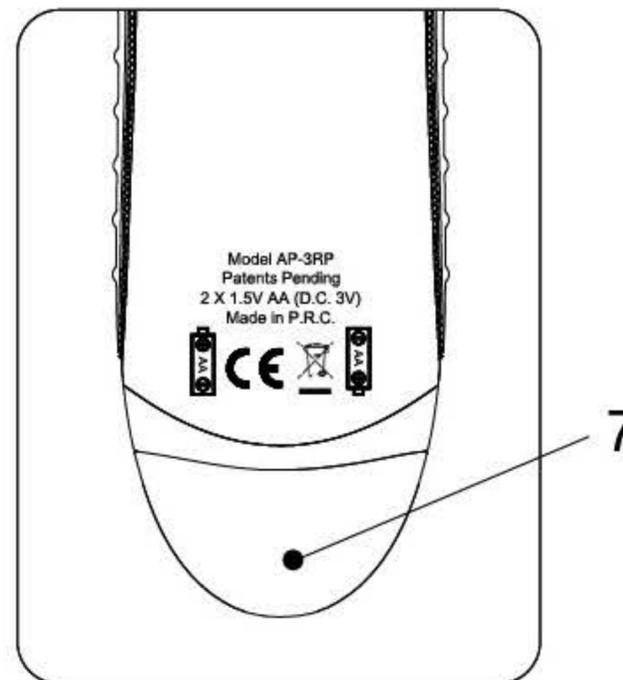

MICRO-Pedi™

Emjoi® Micro-Pedi is an innovative device created to smooth dry, rough and calloused skin on the feet. Unique roller spins 360 degrees, at an amazing 30 times a second, to gently and effectively buff away hard, jagged skin on the soles and heels. **Emjoi® Micro-Pedi** works faster and is safer than metal scrapers, with no blades to cut the skin, and easier to use than other methods, that require constant rubbing to achieve results. It is battery-operated and designed to fit comfortably in your hands. Smooth, silky feet has never been easier.






Each Emjoi® Micro-Pedi roller is recommended for one person only. For hygienic reasons, the Micro-Pedi roller should not be shared with others.

INSTRUCTIONS FOR USE

Before using your Emjoi® Micro-Pedi, please read this instruction manual carefully and completely.

Battery Installation

1. Remove the battery cover (bottom cap of the handle) by turning it 90 degrees to the left until it clicks into place (fig. 1) and then pull it down (fig. 2).
2. Insert 2 AA batteries in the battery compartment following the battery direction sign (fig. 3) as imprinted on the back of the handle.
3. Close the battery cover by turning it 90 degrees to the right (fig. 4) until the cover clicks into place and locks (fig.5).

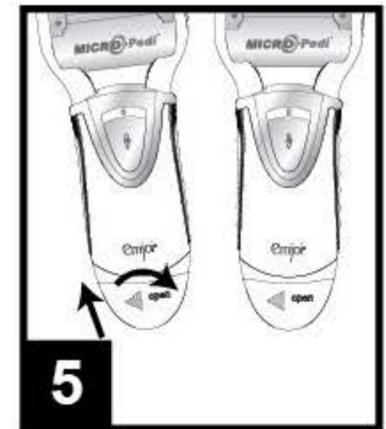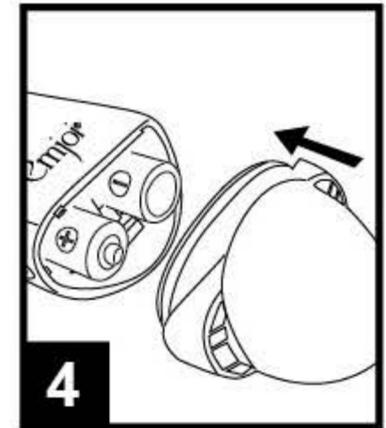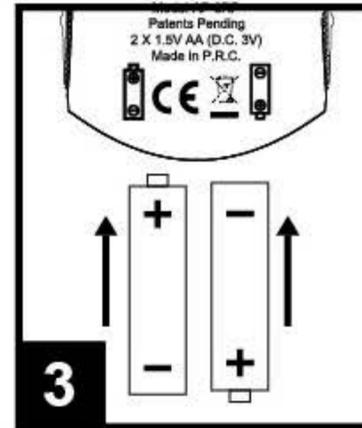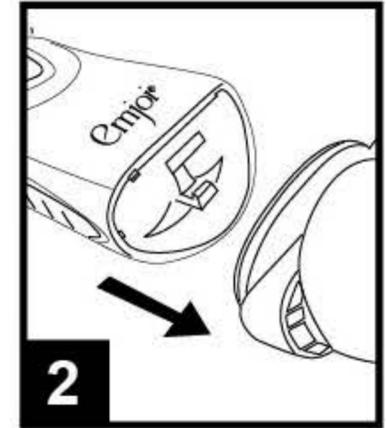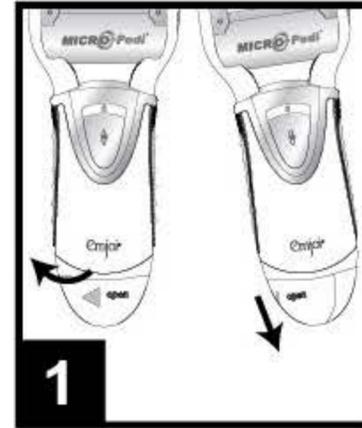

HOW TO USE

Emjoi® Micro-Pedi is designed for smoothing coarse skin on the feet and heels (fig. 6). Do not use this device anywhere else on the body.

Note: Only use Micro-Pedi roller on dry skin.

1. Wash and clean the skin and dry fully.

Note: Do not use on broken or bleeding skin.

2. Make sure the Micro-Pedi roller is placed into the unit. See section of **Changing the Micro-Pedi Roller** on how to remove and replace rollers.
3. Remove the protective cover from the unit and it is ready to use.



4. Turn on the unit by pressing the safety switch lock, in the center of the On/Off switch, and at the same time push the switch up (fig. 7). Gently rotate (back and forth or side to side) over the callous/corn/hard skin area for 2-3 seconds. **Do not press hard- gently glide the unit over desired area.** The unit will stop if it is pressed too hard on the skin. Stop and check if you have achieved the softness you desire. If not, go over the area for 2-3 more seconds and check again. Continue to use this method until the level of smoothness is obtained.

Warning: Stop using if skin becomes sore or inflamed.

Warning: Do not use on skin for more than 2-3 seconds at a time.

5. Turn the unit off (fig 8).
6. Rinse the skin or use a wet towel to get rid of the excess dead particles.
7. To clean the **Emjoi® Micro-Pedi**, please see the Cleaning section.
8. Keep the protective cover on your unit when it is not in use.



CHANGING THE MICRO-PEDI ROLLER

Change the Micro-Pedi roller as needed for best result. Rollers last for a few sessions depending on frequency of use and size of area being treated. Replace as needed.

Replacement Emjoi® Micro-Pedi rollers are available at most retailers that sell the Emjoi® Micro-Pedi unit.

The Micro-Pedi roller contains fragile parts, please handle them with care.



Please take extra care when replacing the Micro-Pedi roller.

Note: Micro-Pedi rollers that are broken or worn could cause damage to the skin.




Replace the **Emjoi® Micro-Pedi** roller when you find any defects to the roller or when you find them less effective.

Make sure the unit is turned off before replacing the Micro-Pedi roller.

1. Pull down the release button, on the side of the unit, while the end cap is facing up and the Micro-Pedi roller will pop out (fig. 9 & 10).
2. After the Micro-Pedi roller pops out, remove it from the cylinder holding the end cap of the Micro-Pedi roller (fig. 11) and dispose of it.

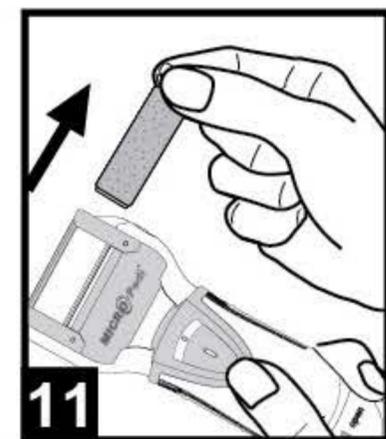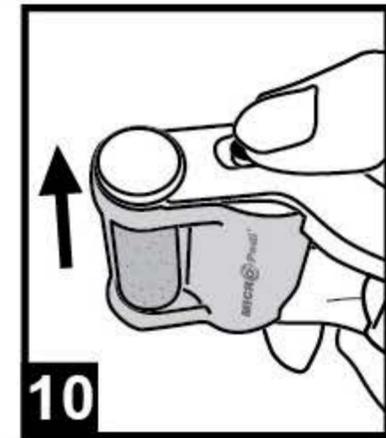

3. Remove the Micro-Pedi roller from the replacement container (fig. 12), (replacement rollers sold separately), by pressing the bottom with your pointer finger and catching it in the palm of your opposite hand.
4. Place a new Micro-Pedi roller into the unit until you hear a 'click' sound. Be sure to insert the new Micro-Pedi roller into the unit with the thinner side going into the unit (fig. 13 & 14).
5. After placing a new Micro-Pedi roller into the unit check that it has been firmly installed with no signs of shaky or loose parts (fig. 15).

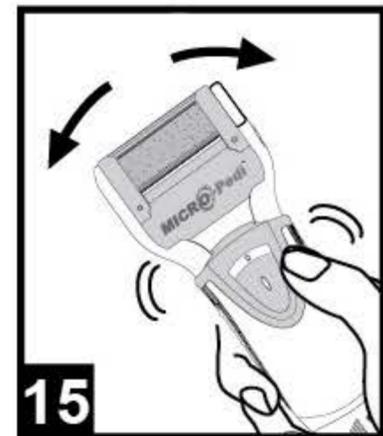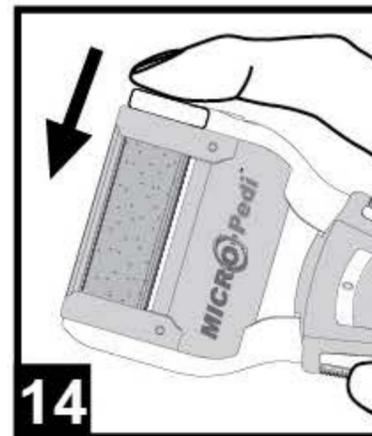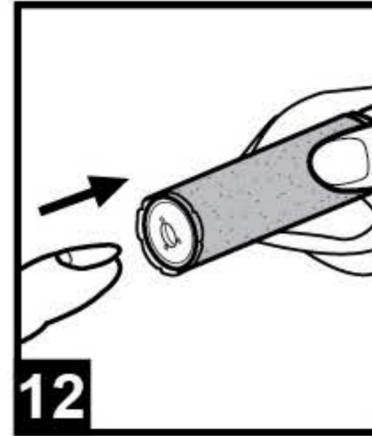

CLEANING

The **Emjoi® Micro-Pedi** is water resistant. It can be rinsed under water for cleaning but it is not water proof. Do not submerge the unit into water on purpose as this may cause defects to the unit. **Micro-Pedi rollers should not be placed into water.**

1. Turn the unit off.
2. Remove the protective cover and the Micro-Pedi roller from the unit.
Note: Micro-Pedi rollers should not be placed into water.
3. Rinse the unit under the tap for 5-8 seconds. Wipe the body of the unit with a clean dry cloth.
4. Replace the protective cover.

WARNING / IMPORTANT

- Before using your **Emjoi® Micro-Pedi** please read all these instructions carefully and completely.
- The **Emjoi® Micro-Pedi** is ONLY for use on the feet and heels. **Do not use this device anywhere else on the body.**
- Use the **Emjoi® Micro-Pedi** only for purposes described in this instruction book.
- Please keep these instructions for future reference.
- This unit is intended for one person. Do not share this unit with others.
- Each **Emjoi® Micro-Pedi** roller is recommended for one person only. For hygienic reason, the Micro-Pedi roller should not be shared with others.
- Do not use if you have diabetes or poor blood circulation.
- If you have any skin conditions or other medical issues, please consult with your physician prior to use.
- Do not use on skin for more than **2-3** seconds at a time.
- Do not use the Micro-Pedi roller when you found worn or nicked edges on the roller.
- Change the Micro-Pedi roller as needed for best result. Rollers last for a few sessions depending on frequency of use and size of area being treated.
- Do not use on skin that is red, irritated, inflamed, infected or blistered.
- Do not use on warts.
- Do not use on sores or open wounds.
- Discontinue use immediately if you experience any pain or irritation.
- Overuse may cause skin irritation.
- Seek medical attention if skin becomes severely irritated or bleeds.
- Only use the **Emjoi® Micro-Pedi** on dry skin.
- This unit is for dry use only.
- Household use only.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the use of the appliance by a person responsible for their safety.

- Do not place or store your unit in direct sunlight. Please store your unit in a cool dry place at a temperature between 0 - 40°C.
- Keep the appliance dry, when it is not in use and after cleaning.
- Do not store your unit where it can fall into a tub or sink of water.
- Never use your unit if it has been damaged or submerged into water.
- Do not use water or any other liquid on the Micro-Pedi rollers.
- Do not turn your unit On if the Micro-Pedi roller is not attached.
- Always turn your unit "Off" when it is not in use and before cleaning it.
- Do not use this unit if the Micro-Pedi roller is damaged as this could lead to injury.
- Do not drop or insert anything into any openings of the unit.
- Do not use your unit outdoors.
- Do not use your unit where aerosol products or oxygen is being used.
- Children should be supervised to ensure that they do not play with this appliance.
- Keep your unit out of reach of children as it is not intended for use by children.
- This appliance is not to be used on, by or near children or invalids.
- To prevent accidents and injury to you or damage to the unit – when your **Emjoi® Micro-Pedi** is in use, keep it away from the hair on your scalp, your eyebrows, eyelashes as well as clothes, brushes, wires, cords, shoelaces, etc.
- If you have any medical conditions, issues, questions or any sensitivity whatsoever we strongly advise you to consult with your health care provider before using this product.
- If you have a pacemaker consult with your health care provider prior to use.
- Do not continually rub the unit on the same area or apply strong pressure.
- Keep the protective cover on your unit when it is not in use.
- Danger of explosion if battery is incorrectly replaced. - Replace only with the same or equivalent type recommended by the manufacturer. Discard used batteries according to the manufacturer's instructions.

Warranty

The **Emjoi® Micro-Pedi** warranty against defects in materials and workmanship lasts one year from the date of purchase.

Keep original receipt for warranty service.

During this warranty period, return the product to Emjoi, Inc. with proof of purchase to FDR Station, PO box 7406, New York, NY 10150. Emjoi will take care of repair, or at their discretion, replacement of a defective unit without charge. For questions regarding this policy you may call toll free 1-888-99-EMJOI (1-888-993-6564) M-F, 9-5 EST.

Important: This warranty does not cover damages resulting from accidents, misuse, alterations, abuse, lack of reasonable care, or performance of services or repairs by an unauthorized service center.

The manufacturers and distributors of the **Emjoi® Micro-Pedi** do not assume responsibility for any damage caused by misuse or use other than that described in these instructions.

Model: AP-3RP

Patents pending

Made in P.R.C.

© Soft Lines International Ltd. 2013

w w w . e m j o i . c o m

PK-IB-EMJ-3RPT-04-0